



KITTY KOLLAR® “E” TUBE SUPPORT

PROTECTS THE WOUND AND ABSORBS DRAINAGE

The placement of esophagostomy “E” feeding tubes is a growing modality in small animal medicine. Many “E” tubes have been successfully in place for years on some patients. The challenge of the “E” tube is its care or maintenance.

The **Kitty Kollar®** was designed as a sensible solution to the everyday patient management of the “E” tube. The adjustable collar fits securely around the animal’s neck. It has a small hole for the “E” tube to pass through. A velcro strap secures or holds the “E” tube flat when not in use. Circular cotton pads with a center hole are provided to be placed between the collar and the skin. They protect the wound and absorb drainage eliminating the need for owners to learn bandage techniques or repeated clinic visits for rebandaging.


FEATURES

- ▶ Machine washable collars provides cleansing of stomasite
- ▶ Easy application and removal
- ▶ Protected against accidental removal
- ▶ Decreases the risk of wound trauma such as scratching
- ▶ Can leave in place for radiographs, chemo, MRI, and other procedures
- ▶ Works for most brands of “E” tubes



SKU	Size
J1085X	X-Small
J1085	Petite
J1085A	Small
J1085B	Medium
J1085C	Large
J1085D	X-Large
J1085K	Small w/E Tube
J1085KS	Petite w/E Tube

All include six protector pads
and laminated instruction card

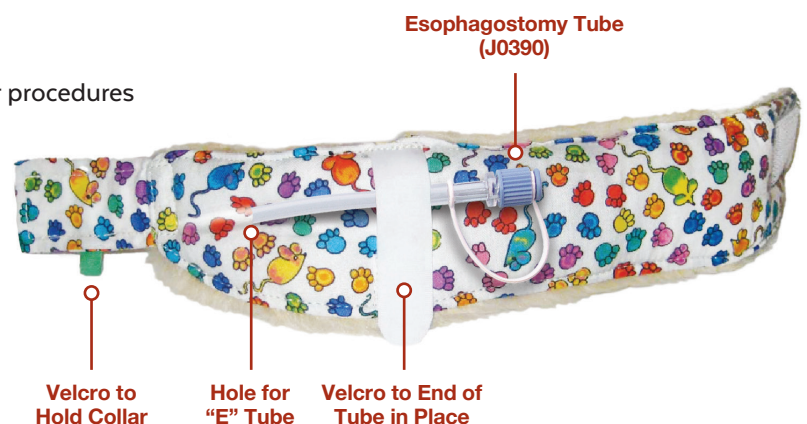

See reverse side for Size Chart, Esophagostomy Tube, and Testimonials



Jorgensen Laboratories, LLC
 800.525.5614 | fax 970.663.5042
www.jorvet.com | info@jorvet.com